

The IonPac NS1 column is an ethylvinylbenzene/divinylbenzene (EVB/DVB) reversed-phase column. This column is ideal for ion-pair or ion suppression separations of hydrophobic ionizable compounds such as organic acids and surfactants. This column can also be used to perform reversed-phase separations. The polymeric nature of the column permits use of acids, bases, and solvents in the eluent.

Column:	IonPac NS1 (10 μ m)	Peaks:	1. Fluoride	5 mg/L
Eluent:	A. 1 mM TBAOH, 1 mM Na ₂ CO ₃ , 10% ACN		2. Chloride	1
	B. 1 mM TBAOH, 1 mM Na ₂ CO ₃ , 20% ACN		3. Nitrite	2
Flow Rate:	1.0 mL/min		4. Bromide	2
Inj. Volume:	50 μ L		5. Nitrate	2
Detection:	Suppressed conductivity, ASRS, external water mode		6. Sulfate	2
Regenerant:	10 mN H ₂ SO ₄		7. Phosphate	3
			8. Thiosulfate	5
			9. Thiocyanate	5
			10. Perchlorate	10

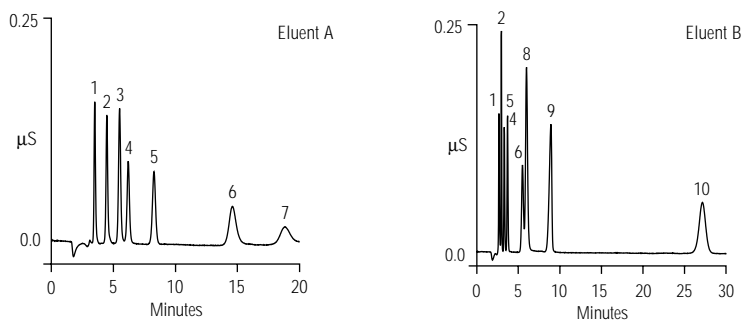


Figure 1 Determination of inorganic anions by ion-pair chromatography.

REVERSED-PHASE OR ION-PAIR CHROMATOGRAPHY

The NS1 column is ideal for Mobile Phase Ion Chromatography (MPIC), which combines ion-pair chromatography with suppressed conductivity detection. This combination can be used to determine high molecular weight ionic analytes that are difficult to separate using ion-exchange chromatography.

VERSATILE DETECTION

Suppressed conductivity detection improves the detection limits for acids or bases by reducing the eluent conductance and maximizing analyte conductance. Alternatively, the NS1 can easily be coupled with UV detection for aromatic molecules.

KEY APPLICATIONS

- High MW aliphatic carboxylic acids
- Anionic and cationic surfactants
- Alkyl sulfates and sulfonates
- Aromatic sulfates and sulfonates
- Quaternary ammonium ions
- Water-soluble vitamins
- Sulfur oxides
- Metal cyanide complexes
- Phenols
- Alkanolamines

WIDE pH RANGE

The IonPac NS1-10 μ m and NS1-5 μ m columns are packed with durable EVB crosslinked with 55 percent DVB. The wide pH range of the IonPac NS1-10 μ m and NS1-5 μ m columns permits the use of eluents from pH 0 to 14. These columns are ideally suited for ion-pair and ion-suppression reversed-phase chromatography—the two techniques that extend reversed-phase chromatography to ionizable compounds. The NS1-10 μ m column is used for routine MPIC separations.

SEPARATE SURFACTANTS AND QUATERNARY AMINES WITH MPIC

MPIC is a good technique for the determination of ionic surfactants and quaternary ammonium compounds in detergents and other matrices. Figures 3 and 4 illustrate the analysis of these compounds using ion-pair chromatography with suppressed conductivity detection.

Aliphatic compounds, which are often nonchromophoric, can easily be determined using MPIC. When combined with ion-pair chromatography, suppressed conductivity is a sensitive, universal detection mode for analytes such as aliphatic carboxylic acids, sulfonic acids, and quaternary ammonium ions. For aromatic acids or bases, UV detection provides optimum sensitivity.

As illustrated in Figure 5, MPIC can be used by food and cosmetics manufacturers for the determination of additives such as alkanolamines. In the pharmaceutical industry, the technique can be used for a variety of non-chromophoric, ionic species.

USE THE NS1 COLUMN FOR ION SUPPRESSION

Ion suppression is a technique used to suppress the ionization of weak acids to enhance reversed-phase retention. As illustrated in figures 6 and 7, ion suppression can be used to separate a variety of compounds, including aliphatic carboxylic acids and substituted benzoic acids.

CHEMICAL SUPPRESSION WITH CONDUCTIVITY DETECTION

Anion-MPIC

The eluent suppressor used for anion-MPIC is the Dionex ASRS, which is designed to suppress ammonium or tetraalkylammonium ion-pair reagents as well as to suppress conventional anion-exchange eluents. The eluent suppressor used for cation-MPIC is the Dionex CSRS, which can also be used to suppress conventional cation-exchange eluents.

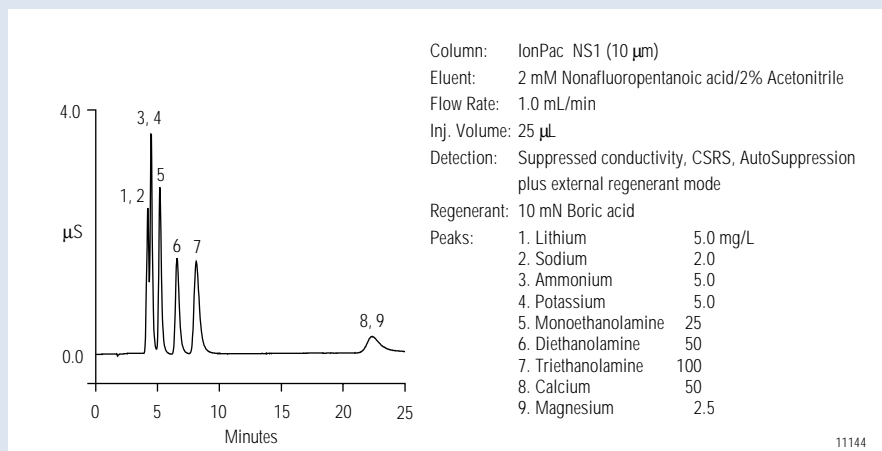


Figure 5 Separation of alkanolamines on the NS1 column using ion pairing with suppressed conductivity detection.

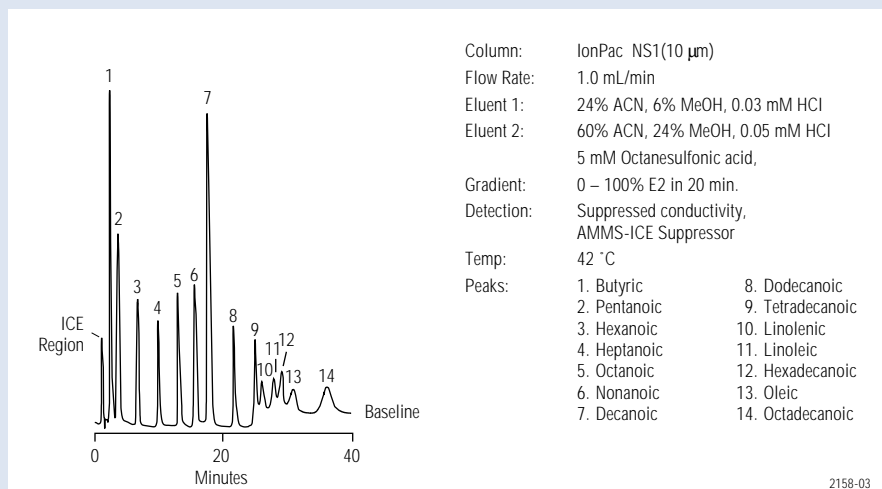


Figure 6 Determination of aliphatic carboxylic acids using the IonPac NS1 column.

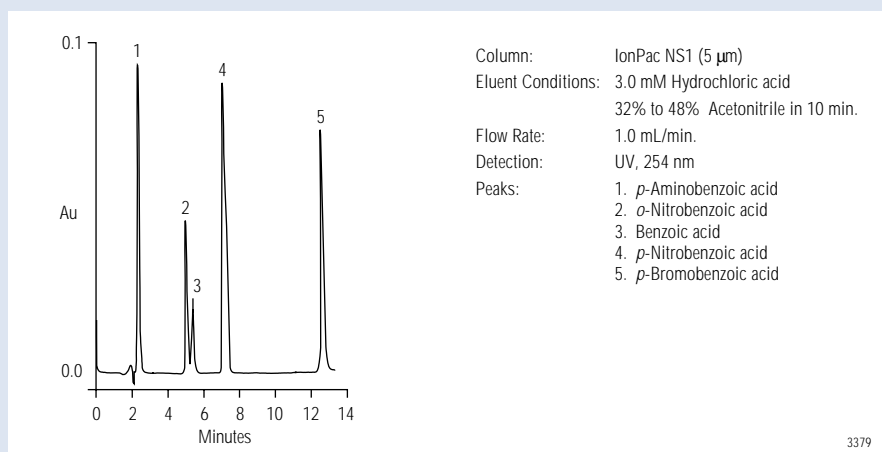


Figure 7 Determination of substituted benzoic acids using ion-suppression chromatography.

Using the ASRS, the water regenerant undergoes electrolysis to form hydronium ions.

The hydronium ions neutralize the basic eluent, while ion-pair counterions of the basic eluent or sample move across the membrane to waste. The result of this suppression process is a sample of anions paired with highly conductive hydronium counterions resulting in optimum conductance in a water eluent.

In anion-MPIC, the quaternary ammonium counterions have a very high affinity for the cation-exchange membrane. The sulfuric acid regenerant is added to effectively displace the quaternary ammonium ion-pair reagent from the membrane, lowering the resistance in the suppressor.

Cation-MPIC

The Cation Self-Regenerating Suppressor (CSRS) operation is very similar to the ASRS. The CSRS is used for the neutralization of acidic eluents containing anionic ion-pair reagents.

Ion Suppression

The AMMS-ICE Suppressor is used for the acidic eluents used in ion suppression chromatography of organic acids. Potassium hydroxide is typically used as the regenerant for the AMMS-ICE suppressor. The AMMS-ICE contains cation-exchange membranes, which exchange the hydronium ion of the acids for potassium ion, thus causing complete ionization of the acids and allowing conductivity detection of the analytes.

SPECIFICATIONS

Analytical Column Dimensions

IonPac NS1-10 μm , 4 x 250 mm

IonPac NS1-5 μm , 4 x 150 mm

Guard Column Dimensions

IonPac NG1-10 μm , 4 x 35 mm

Maximum Operating Pressure:

<27.6 MPa (<4000 psi)

Mobile Phase Compatibility:

pH 0-14; 1-100% HPLC solvents

Column Construction:

PEEK with internal style 10-32 threaded ferrule-style end fittings. All components are nonmetallic.

Substrate Characteristics:

Particle Diameter: 5 and 10 μm

Crosslinking (%DVB): 55%

Hydrophobicity: High

ORDERING INFORMATION

In the US, call 1-800-346-6390 or contact the Dionex Regional Office nearest you. Outside the US, order through your local Dionex Office or Distributor. Refer to the part numbers below.

Analytical and Trap Columns

IonPac NS1-10 μm Analytical Column (for routine MPIC determinations) (4 x 250 mm) P/N 35321

IonPac NG1-10 μm Guard Column (4 x 35 mm) P/N 39567

IonPac NS1-5 μm Analytical Column (for high resolution MPIC determinations) (4 x 150 mm) P/N 39568

IonPac ATC-1 Anion Trap Column (for anion gradient operation) (9 x 24 mm) P/N 37151

IonPac CTC-1 Cation Trap Column (for cation gradient operation) (9 x 24 mm) P/N 40192

Suppressors

ASRS (4-mm) Anion Self-Regenerating Suppressor for anion ion pairing P/N 43189

CSRS (4-mm) Cation Self-Regenerating Suppressor for cation ion pairing P/N 43190

AMMS-ICE Micromembrane Suppressor for anion ion suppression P/N 37107

MMS/SRS External Regenerant Kit P/N 38018

Ion Pairing Reagents

Tetrabutylammonium Hydroxide, 500 mL 0.1 M TBAOH (MPIC-AR1) P/N 35360

Tetrapropylammonium Hydroxide, 500 mL 0.1 M TPAOH (MPIC-AR2) P/N 35363


Hexanesulfonic Acid, 500 mL 0.1 M HSA (MPIC-CR1) P/N 35361

Octanesulfonic Acid, 500 mL 0.1 M OSA (MPIC-CR2) P/N 35362

Regenerant Reagents

Anion Regenerant Solution, for MPIC, 4 x 200 mL 0.5 N H₂SO₄ 20X concentrate P/N 37164

IonPac is a registered trademark and AutoSuppression is a trademark of Dionex Corporation.

 Printed on recycled and recyclable paper with soy-based ink.

Dionex Corporation
1228 Titan Way
P.O. Box 3603
Sunnyvale, CA
94088-3603
(408) 737-0700

Dionex Corporation
Salt Lake City Technical Center
1515 West 2200 South, Suite A
Salt Lake City, UT
84119-1484
(801) 972-9292

Dionex U.S. Regional Offices
Sunnyvale, CA (408) 737-8522
Westmont, IL (630) 789-3660
Houston, TX (713) 847-5652
Smyrna, GA (770) 432-8100
Marlton, NJ (609) 596-0600

Dionex International Subsidiaries
Belgium (015) 203800 Canada (905) 844-9650 France (1) 39 46 08 40 Germany (06126) 991-0 Italy (6) 66030052
Japan (06) 885-1213 The Netherlands (076) 57 14 800 Switzerland (062) 205 99 66 United Kingdom (01276) 691722
* Designed, developed, and manufactured under an NSAI registered ISO 9001 Quality System.
<http://www.dionex.com>

